




Accuris Mini Balances are extraordinarily compact, easy to use, and convenient for a variety of laboratory weighing applications. With a footprint of only 2.6 x 3.9 inches, these Mini Balances offer exceptional accuracy and linearity through their full weighing ranges.


Available with either a 100g or 500g capacity, these digital balances feature an easy to read, backlit, liquid crystal display. Buttons for ON/OFF, MODE and TARE are located on the front for easy access. The software includes a calibration function for easy in-lab calibration with a 100g or 500g calibration weight. Calibration weights are available separately.

Designed for portability, the Accuris Mini Balances are powered by 2 x AAA batteries (included). After one minute sitting idle, the auto shut off feature turns off power to conserve battery life. Ideal for lab or field use, each balance includes a snap on, protective cover with instructions conveniently printed on the underside

- Easy-to-read backlit liquid crystal display
- Calibration function
- Includes 2 x AAA batteries
- Snap on, protective cover
- Compact, less than 3 in. wide

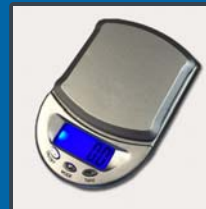

Mini 500g

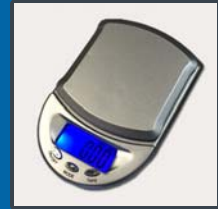

Mini 100g



Protective cover (included)

Model	W4000-100	W4000-500
Max Capacity	100 grams	500 grams
Readability	0.01 g	0.1 g
Repeatability	0.01 g	0.1 g
Linearity (+/-)	0.01 g	0.1 g
Calibration	External	
Pan Size	2.1 x 2.25 in / 5.3 x 5.7 cm	
Weighing units	g, ct, dwt, gn	
Display	backlit LCD	
Power	2 x AAA batteries	
Exterior Dims.	2.6 x 3.9 x 0.75 in / 6.6 x 9.9 x 1.9 cm	

Accuris Instruments offer a full line of Compact, Precision and Analytical Balances.



Visit: www.accuris-usa.com

Ordering Information

Item No.	Description:
W4000-100	Accuris™ Mini Balance, capacity 100 grams / 0.01g readability
W4000-500	Accuris™ Mini Balance, capacity 500 grams / 0.1g readability
W1000-100	Accuris™ Calibration Wight, F1 Grade, 100g
W1000-500	Accuris™ Calibration Wight, F1 Grade, 500g